



Disposable Injection Needles

The widest choice of both U.G.I and L.G.I needles available in the UK, designed to ensure safe, effective injection with the greatest ease of operation.



Diagmed Healthcare Ltd is a specialist distributor of world leading medical technology.

Diagmed's mission is to capitalise on our experience and assets by identifying and representing unique and innovative product lines that benefit patients and medical staff.

We fully understand exactly how to distribute medical/surgical devices within the UK health care markets. We share a good working relationship with a very large customer base and fully support our customers before and after the point of sale.

Our team of product specialists cover the whole of the UK, undergo regular training on product knowledge and have the full technical support from a team of administrators at head office.

Sclerotherapy Needles




Sure Stop

The unique design of the Diamed Sclerotherapy Needle ensures safe, high quality endoscopic procedure with the greatest ease of operation available. Design features include:

The Sure Stop Handle

- Permits one handed extraction of the entire instrument from the scope
- Sure Stop lock prevents separation of inner and outer sheath during retraction
- Protects equipment from perforation: built in advance and retraction stops provide maximum needle extension and prevent over-withdrawal of needle into outer sheath
- Eliminates guesswork about the proper position of the retracted needle
- Protects patient by eliminating the possibility of accidental penetration during multiple insertions

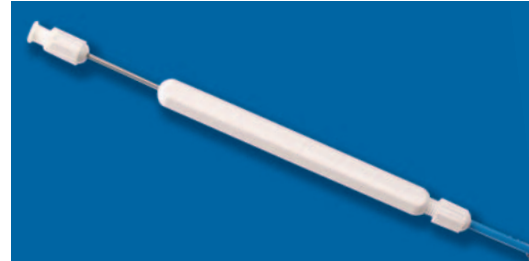

The Sure Stop Sheath

- Assures consistent 5mm needle extension beyond the distal tip, preventing over penetration.
- 1.8mm and 2.3mm diameter sheaths for a wider range of applications
- Transparent sheaths permit visual identification of sclerosant flow, injectant or blood

The Sure Stop Handle

- Standard 200cm working length covers every procedure and reduces the need for different length inventories
- Available in three sizes, 25 gauge needle for a small, precise puncture and larger 23 gauge and 21 gauge bores for higher flow rate and faster injection, allowing bulkier injectant to be used.

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
DIS2002	23	1.8	5	200
DIS2004	23	2.3	5	200
DIS2006	25	1.8	5	200
DIS2008	25	2.3	5	200
DIS2010	21	2.3	5	200
DIS2014	25	1.8	5	230
DIS2015	25	2.3	5	240
DIS2016	23	2.3	5	240
DIS2018	25	2.3	5	160

Optima Sclerotherapy Needle

One of the superior ranges of needles on offer, this needle features a locking cap to ensure that the needle is clearly either deployed or retained to avoid sticking and subsequent damage either to scopes or to patients. Once locked the needle will remain deployed until retraction without user effort.



- Smooth distal tip hub – helps prevent sheath perforation
- Enhances maneuverability
- Clear outer teflon sheath enables user to monitor needle and movement of inner sheath
- Maximises needle projection in retroflex positions
- Eliminates kinking

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
DIS2002LC	23	1.8	5	200
DIS2004LC	23	2.4	5	200
DIS2008LC	25	2.4	5	240
DIS2010LC	21	2.4	5	200
DIS2016LC	23	2.4	5	240
DIS2017LC	23	2.4	5	350
DIS2026LC	23	2.4	4	200

Optima Plus Sclerotherapy Needle

Our superior needle with the added benefit of a metal sheath for greater flexibility especially for use in ERCP. Economically priced against bigger brand names.

- Teflon coated stainless steel outer sheath
- Metal inner jacket
- Smooth ball top supporting needle during articulation
- Enhances maneuverability
- Maximises needle projection in retroflex positions
- Eliminates kinking
- Locking handle (not spring locked)

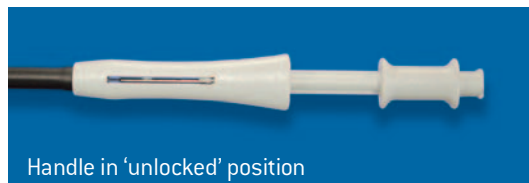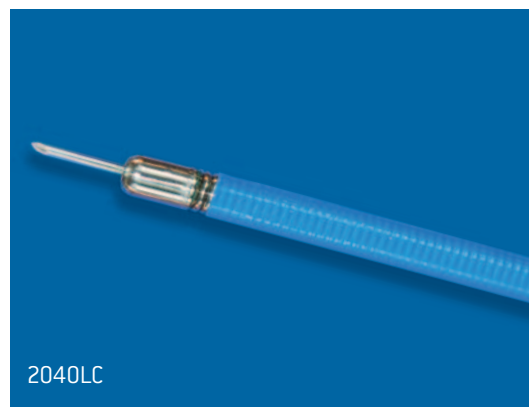


Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
DIS2040LC	25	2.4	5	230

Economy Range of Injection Needles

Sure Click

- Click Lock System
- Distal Guide Ring
- Handle permits one handed usage

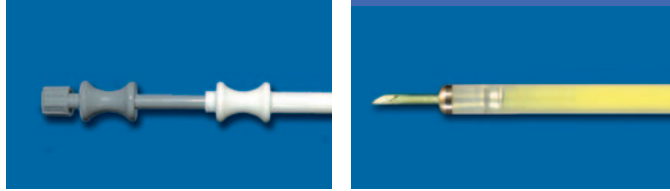

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
09-10-72-32-12AS	23	2.35	5	230
09-11-02-32-12AS	19	2.35	5	230

Eco-Click Injection Needle

- Blue teflon inner and outer tube
- Stainless steel needle
- 4,5,6 and 8mm needle lengths available
- Available with 19/21/23/25 needle gauge
- Available for Gastric and Colonic procedures

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
AF-120-25-4	25	1.8	4	120
AF-165-25-5	25	2.4	5	165
AF-230-19-5	19	2.4	5	230
AF-230-21-5	21	2.4	5	230
AF-230-25-5	25	2.4	5	230

A Guide to Product Codes and Needle Dimensions

Example Product Code AF-XXX-YY-Z

XXX - Working Length (cm) YY - Needle Gauge Z - Needle Length (mm)

All of these needles are priced for economy but designed and manufactured for performance.

Single Use Injection Needle

Diagmed's single use injection needle is cost effective for a range of mucosal therapy procedures.

Design features include;

- Click lock handle
- Safe and secure deployment and retraction of needle

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
SN-1623-22	22	2.3	5	230

Injection Needle for Hystoacryl (Glue Injection)

Diagmed's injection needle for Hystoacryl offers the following features;

- Click Lock Handle for secure deployment
- Luer Lock Syringe for secure attachment
- Metal Tipped Distal End for greater patient, scope and user safety
- Wide bore needle for greater ease of injection



Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
99-09-30-72-32-07	22	2.3	5	230

Medi Jet Sclerotherapy Needles

Needle for effective sclerotherapy of varices within the GI tract by means of vascular or submucosal injection - in case of acute bleedings or as a preventive treatment.

Design features include;

- Injection needle with locking mechanism
- Catheter tip with distal tapering (filiform)
- Transparent catheter to aid visualisation of injectant.

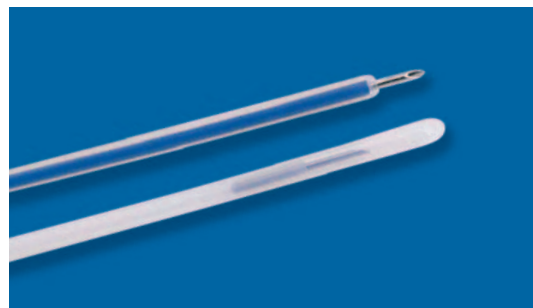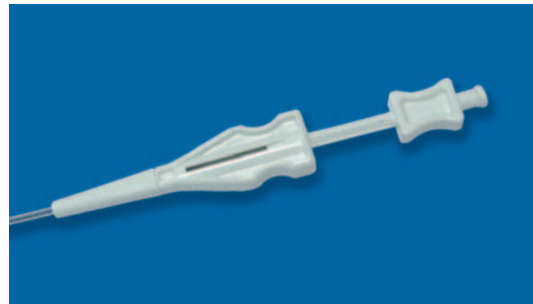

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
E707230	22	2.3	5	230
E710180	20	2.3	5	180
E710230	20	2.3	5	230

Transbronchial Needle

Diagmed's Transbronchial Needle ensures safe, high quality endoscopic procedure with the greatest ease of operation available. Design features include:

The Sure Stop Handle

- As per features of the Sclerotherapy Needle range

The Sure Stop Sheath

- Assures reproducible 13mm needle extension beyond the outer sheath, preventing over penetration
- 1.8mm outer catheter diameter fits most bronchoscopes and allows full scope augulation
- The 150cm catheter length will allow complete bronchoscopic diagnosis
- Eliminates the need for a stiffening stylet

The Sure Stop Needle

- Stainless steel 21 gauge needle is designed for maximum specimen retrieval with a very small puncture
- Needle aspiration is particularly effective in diagnosing malignancy in patients where traditional sampling instruments cannot fully penetrate the mass

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
DIS2102	21	1.8	13	150

Transbronchial Aspiration Needle

The Transbronchial Aspiration Needle is used in the bronchoscopy procedure to take tissue samples from lymph nodes. The process of taking a lymph sample requires the Doctor to undertake a 'blind' procedure where he punctures the 'narrowing' bronchial canal with the transbronchial aspiration needle, then pushes the needle into the lymph node.

A syringe (usually 20ml) is attached and aspiration is then performed to suck tissue sample up into the needle. A needle with a side port is ideal as this generally gives the Doctor a better tissue yield.

Product Code	Needle Gauge	Sheath Diameter (mm)	Needle Length (mm)	Working Length (cm)
ASP-15-18	25	1.8	15	120

