



Foreign Body Retrieval Forceps



Diagmed offers an extensive choice of disposable foreign body graspers designed to handle any clinical situation encountered during endoscopy procedures.

Disposable Foreign Body Removal Forceps

Diagmed Healthcare retrieval forceps have several unique advantages - the ability to grasp objects of up to 3cm in size, and are made from stainless steel. This ensures durability and precision performance. Our forceps are manufactured using the same high quality materials as reusable forceps.

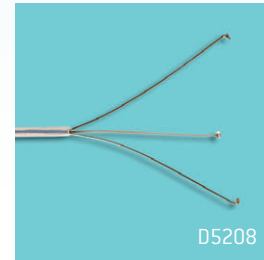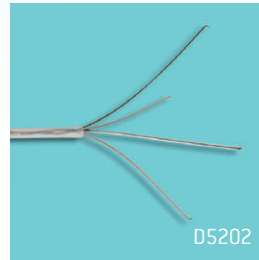


Product Code	Description	Sheath Diameter	Working Length	Open Diameter
D5202	Four Prong Grasping Forcep	2.3mm	240cm	2.5cm
D5204	Four Prong Grasping Forcep	1.8mm	180cm	2.5cm
D5208	Three Prong Grasping Forcep	2.3mm	240cm	2.5cm

Disposable Rat Tooth and Alligator Jaw Forcep

Both the Alligator and Rat Toothed Forceps are used to handle difficult to hold objects. Their design of interlocking teeth and serrated edges enables a greater grip on objects and less slippage reducing repeated grasping.

- Teflon coated for easy passage through scope
- Excellent flexibility
- Strong grip for plastic stent removal

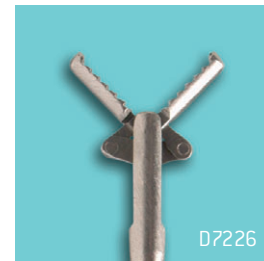


Product Code	Description	Sheath Diameter	Working Length	Jaw Length
D7224	Rat Tooth Jaw	2.4mm	230cm	3.5mm
D7226	Rat Tooth Alligator Jaw	2.4mm	230cm	8.0mm

Disposable Baskets

Mini Catch


Originally designed for ERCP procedures the Mini Catch can also be used for Polyp/Foreign Body Retrieval due to its length and flexible and easy to use basket size and shape.

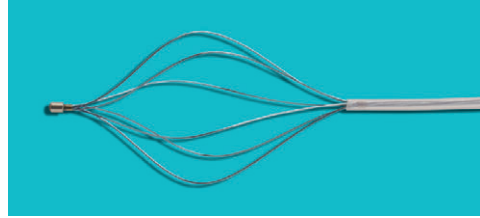

Product Code	Cannula Diameter (mm)	Working Length (cm)	Basket Size (mm)
99-03-12-12-09	1.0	230	15/30

Retrieval Basket

Diagmed Healthcare's four wire basket is constructed with high strength stainless steel wire for maximum strength even when fully extended. The disposable retrieval basket may be used to facilitate the retrieval of foreign bodies, stones or polyps. The basket incorporates four wires, which are extremely flexible allowing the endoscopist to force the basket over the object to be removed, yet maintaining a small cage for retention of even the smallest object.

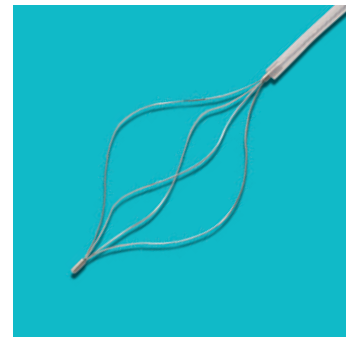

Product Code	Description	Sheath Diameter (mm)	Working Length (cm)	Basket Dimensions
D5210	Disposable Retrieval Basket	2.3	240	20mm wide x 4cm length