


Infinity™ sampling device



Looking to collect quality samples during ERCP?

The Infinity™ sampling device is purpose built for collecting substantial and quality samples from strictures in the biliary duct.

The Infinity™ device features...

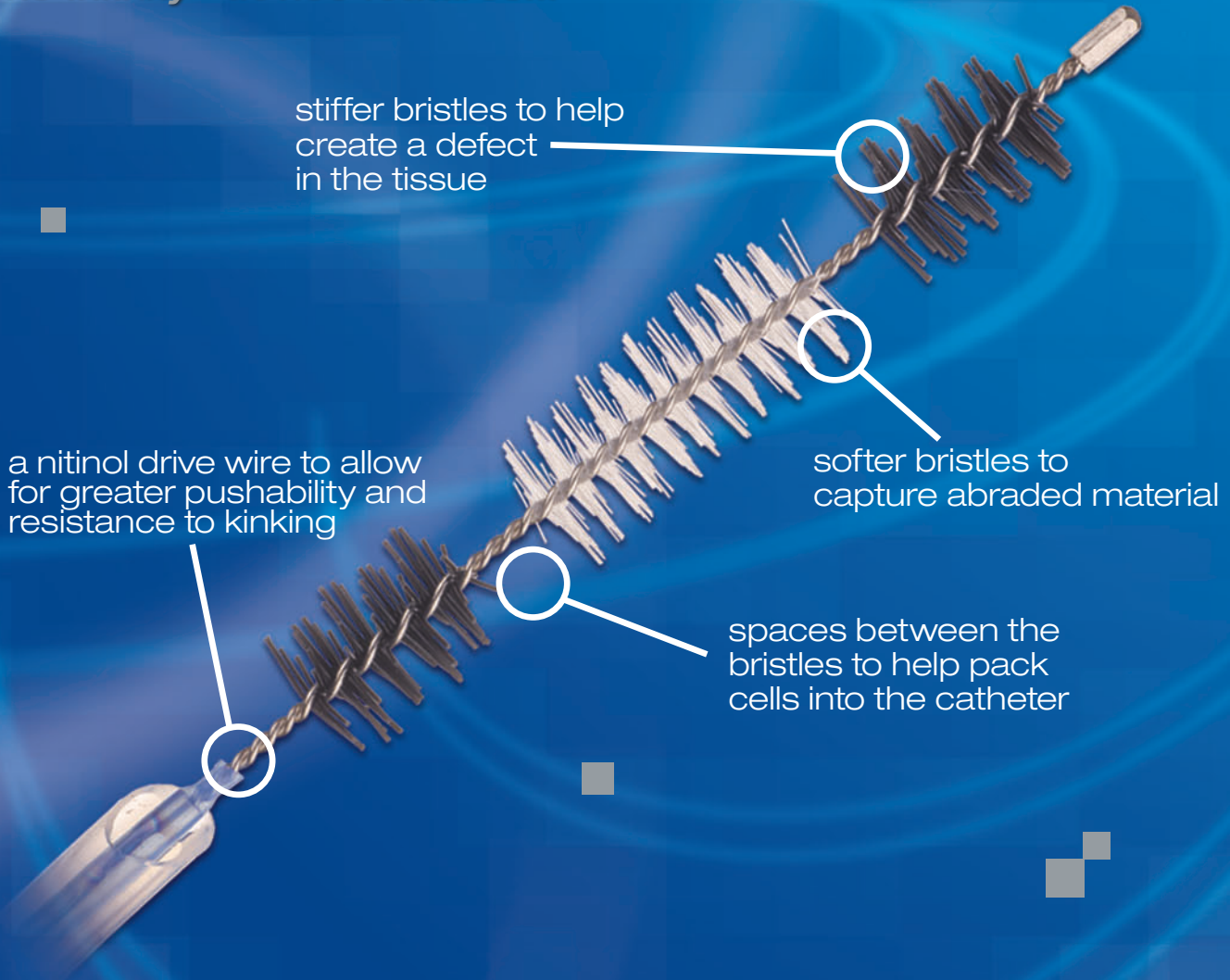

stiffer bristles to help create a defect in the tissue

a nitinol drive wire to allow for greater pushability and resistance to kinking

softer bristles to capture abraded material

spaces between the bristles to help pack cells into the catheter

For more information on the Infinity™ sampling device, contact US Endoscopy today.

product number	sheath diameter (mm)	length (cm)	sterile	units/box
00711652	3 / 9FR	200	yes	5

Distributed by:



listening...and delivering solutions®

5976 Heisley Road
Mentor, OH 44060

phone 440 / 639.4494

fax 440 / 639.4495

customer service 800 / 769.8226

www.usendoscopy.com

760746 Rev.A